

RAVEN[®] LATERAL

Lumbar Lateral Plate System



RAVEN[®] LATERAL

Lumbar Lateral Plate System


The ChoiceSpine Raven[®] Lateral Lumbar Plate System is a two-screw plate for lateral lumbar spinal fusion. Raven Lateral is a two-screw plate with an integrated cam lock mechanism for screw blocking. The Raven Lateral Plate System has simple instrumentation that is compatible with the ChoiceSpine VEO[®] Lateral Access System.

Features:


- Compatibility with our VEO Lateral Access System
- Integrated cam locking mechanism that provides a visual, tactile and audible confirmation of screw blocking.
- 10 degrees of screw angulation
- Plate Heights: 15mm – 23mm in 2mm increments
- Screw Lengths: 35mm – 60mm in 5mm increments
- Screw diameters: \varnothing 5.0mm & \varnothing 5.5mm
- Streamlined instrumentation



Compatible with
VEO[®] Lateral
Access System

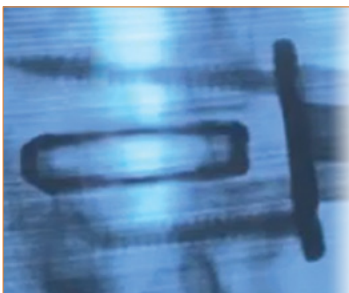

Raven Lateral with
TigerShark[®] Lateral Interbody

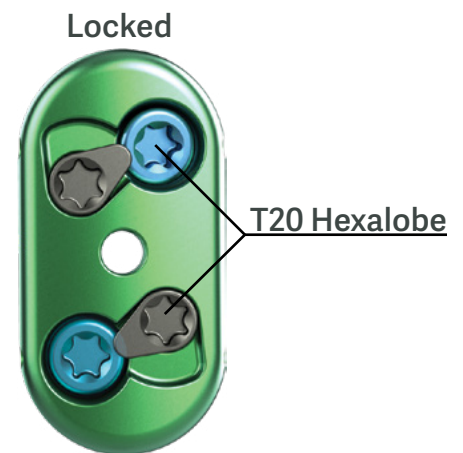

400 Erin Drive, Knoxville, TN 37919 | O: 865.246.3333 | F: 865.246.3334 | choicespine.com